

Invented for life



BOSCH

Vivasuite for Vivatmo

The digital ecosystem for your Vivatmo *pro* device

✓ Vivatmo from Bosch

The intelligently connected FeNO measuring system.



✓ **ivatmo**

What is Vivasuite?

Vivasuite is the digital healthcare ecosystem by Bosch. It is not a laboratory information system, but a cloud based device management tool, that enables you to manage your devices at different locations and receive system updates to ensure availability of the most recent features. Vivasuite works independently from your LIS/HIS system. Vivasuite runs in the Bosch managed Cloud and applies highest standards regarding IT security and data privacy.



- ▶ **Get the newest features**
- ▶ **Minimize system downtime**
- ▶ **Track your device performance**
- ▶ **Export data for statistics**
- ▶ **Control your data**
- ▶ **Manage multiple devices at remote settings**



Connect your device with Vivasuite

Manage software updates

- Stay updated on new device features
- Track your software history


Technical data

- Supply Vivasuite with technical information to ensure continuous device performance

Data security





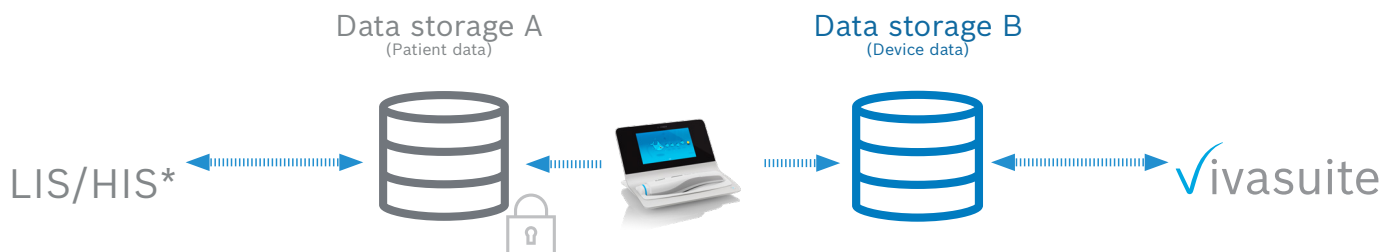





- No access to patient data
- No remote access to the device
- Separate data storages for different purposes
- Vivatmo exchanges solely technical information with Vivasuite

Vivatmo data storage separation



*Vivatmo pro ready for third party plugins via GDT and HL7-Interface

Advantages of Vivasuite

 Receive remote software updates	 Know when your devices were last synchronized	 View all devices in one dashboard	 View device information and general device data	 Check for remaining measurement trials
 Bosch managed Cloud device	 Mobile device monitoring	 Vivasuite user management	 Save time with less "on-device" work	 Get statistics relating to your device usage

Bosch Healthcare Solutions GmbH Stuttgarter Str. 130, Waiblingen, Germany

Vivasuite is not classified as a medical device. Subject to technical modifications.
Not all products are available in all regions. Ask your local sales representative for availability in specific markets.